



Be READY

Casualty Care Rescue Randy

Nasco
HEALTHCARE

Ready to treat. Ready to act.

Introducing **Casualty Care Rescue Randy**

Full-body manikin that is uniquely designed for enhanced realistic training on the procedures that treat the 3 most preventable causes of death:

1. Massive Bleeding
2. Tension Pneumothorax
3. Airway Obstruction

- Compliant with TCCC/TECC scenarios
- Rugged, reliable, and hyper-realistic
- 100% hand-made in the USA



Casualty Care Rescue Randy: Life-like training solutions for real-life situations.

Massive Hemorrhage


- Femoral and Brachial extremity wounds
- Inguinal and Carotid junctional wounds
- Treatment Interventions
Extremity and junctional tourniquets, wound packing and pressure bandage

Airway Management

- Nasopharyngeal (Nose)
- Oropharyngeal (Mouth)
- Cricothyroidotomy (Throat)

Tension Pneumothorax

- Needle Decompression at the 4th/5th intercostal space in the anterior axillary line
 - Needle Decompression at the 2nd/3rd Intercostal space in the mid-clavicular line
- ### Penetrating Chest Wound
- Allows for occlusive dressing

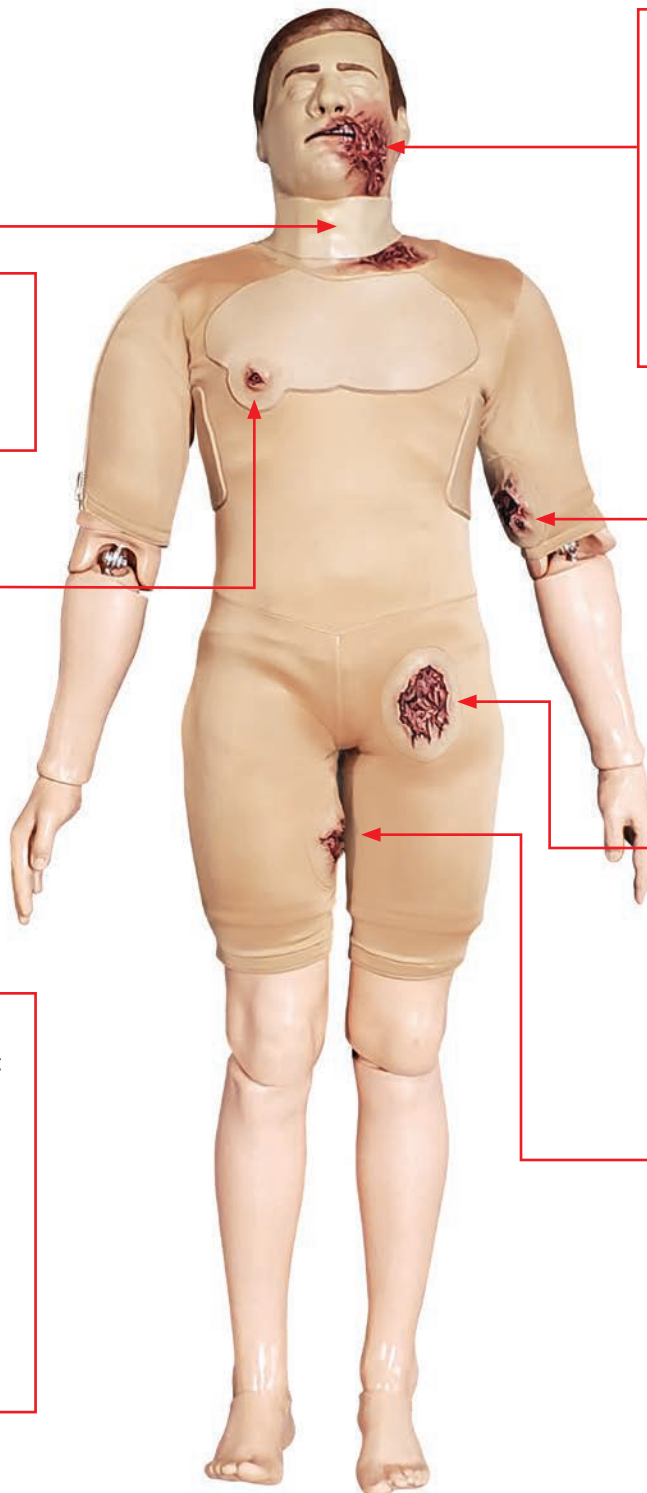

Casualty Care Rescue Randy product profile

Product Name	Casualty Care Rescue Randy - powered by Strategic Operations Hyper-Realistic® technology
Addresses Training Needs/Curricula	TCCC - Tactical Combat Casualty Care for all service members TECC - Tactical Emergency Casualty Care for all EMS/EMT/LEO professionals
Configuration	<ul style="list-style-type: none"> • Full-body manikin with articulated joints • Weight-distributed to mimic the human anatomy • Holds 3-4 liters of blood • Pumping system with regulator simulating human blood flow pressure (2-3) psi • Manual foot operated blood pumping system allows 3 wounds to bleed simultaneously • Optional Electric Deluxe blood pumping system allowing 4 wounds to bleed simultaneously • Comes with blood powder (makes 5 gallons) • Three (3) repairable neck skins • Ten (10) repairable tracheas • One (1) repair kit
Dimensions	Male at 72" tall x 25" wide x 12" deep (165.1 cm x 63.5 cm x 30.48 cm)
Product Ship Weight	122 lbs. (2 boxes)
Product Code/SKU	149-6075



Nasco Healthcare Inc.
 16 Simulaids Drive
 Saugerties, NY, 12477 USA
 Phone: 1.800.431.4310
 Email: info@nascohealthcare.com
www.nascohealthcare.com